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P.I.R.E: <=0.49 dBm.



TERRAILLON CONNECTED BODY COMPOSITION ANALYZER

Thank you very much for selecting the Terraillon Glass Body Fat Analyzer BODY COLOR CONNECT.

Don't use this product if you wear a pacemaker or implantable electronic device • Any information provided by this device is in no way meant to treat, cure or prevent any disease or illness from happening.

PRODUCT OVERVIEW

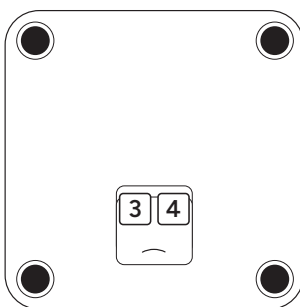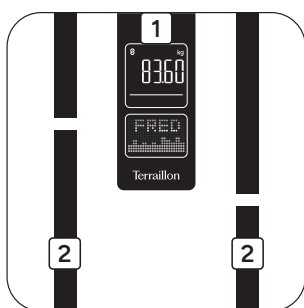- 1** Big colour display
- 2** Electrodes
- 3** Battery compartment
- 4** Unit button
(under the battery compartment cover)

FIRST USE

- Insert the batteries (4xAAA).
- You will need your smartphone to install your bathroom scales.
- Download the Wellness Coach - MyHealth app, then open it.

- Follow the instructions that appear on your smartphone. Ensure you enter all the information requested.
- Your smartphone must remain near to the bathroom scales throughout the installation process.
- You can connect to your scale as the main user or as a guest. Please note: only the main user of the scale will be able to change its settings and delete or add guests.

If you have any difficulties or need more information, you can find our help centre by scanning this QR code with your smartphone. If you cannot find the answer to your question, you can also contact our After Sales Service: support.uk@terraillon.com


TERRAILLON CONNECTED BODY COMPOSITION ANALYZER

WEIGHING OPTIONS

- **Pregnancy mode:** when pregnancy mode is selected, body mass is not taken into account. Only your weight will be measured.

After the first weighing session, to change your options at any time, simply go to the Device tab of your MyHealth application, select , and then select or deselect the option you want to change.

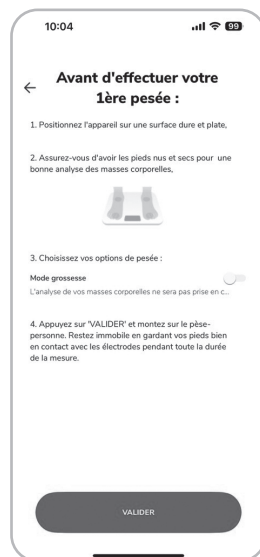START MEASURING

1. You must have bare feet to be able to analyse your body mass.
2. Open the MyHealth application.
3. Step onto your bathroom scale. Maintain contact with the electrodes and wait for the measurement to complete. The data is transmitted directly to the app.
4. If the application is not open, the data will not be saved on the bathroom scale.
5. In the case of weighing without the app open, and if the user is recognized (+/- 2Kg), the scale will display weight, body fat percentage and BMI twice in succession before switching off.

Be sure to always weigh yourself on a hard and flat surface.


MULTI-USER MANAGEMENT

When the device is paired, your MyHealth application will ask you to become the main user. There can only be one main user of the device. This user has the rights to modify certain settings, such as the scale configuration.

To invite a new user while remaining the main user of your BODY COLOR CONNECT, you can follow the steps below to display the QR code needed to pair up the invited user.

The guest user can add a device by scanning the QR code (add device / Scale / Connection via Scan code). The guest user's profile is then associated with the bathroom scale.

TERRAILLON CONNECTED BODY COMPOSITION ANALYZER

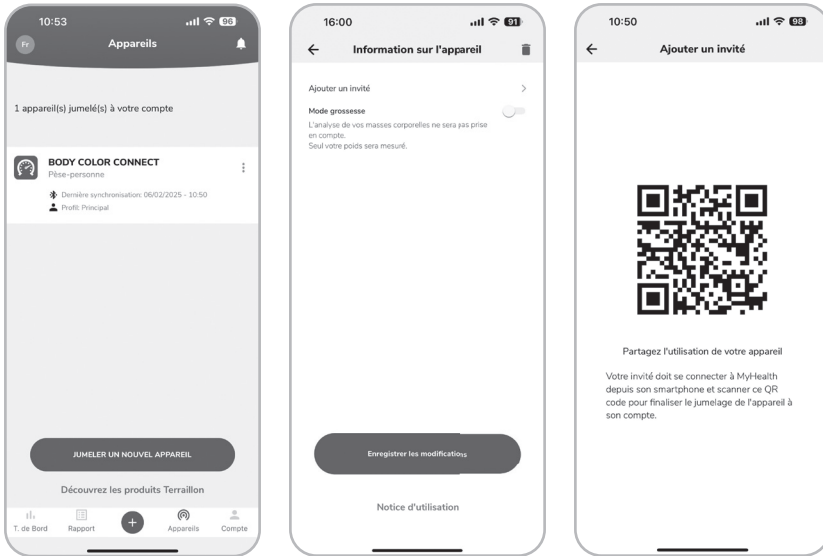MAINTENANCE

Use water only. DON'T use propellant or other chemicals.

WARRANTY & ENVIRONMENT PROTECTION

This product is guaranteed for 2 years against material and manufacturing defects. During this period, such defects will be repaired free of charge (proof of purchase must be presented if the scale is under warranty). This warranty does not cover damage resulting from accidents, improper use or negligence. If you have a complaint please contact the store where you purchased the scale first.

Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice.

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iPhone is registered trademark of Apple Inc.
Android is a trademark of Google Inc.
The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.
Other trademarks and trade names are those of their respective owners.
Terraillon SAS France & Headquarters
Terraillon Asia Pacific Ltd